4/27/22, 9:11 AM PCR preperation

PCR preperation

2022-04-27

Dict{Any, Any} with 7 entries:

"FileName" => "PCR.fasta"

"Cycles" => 30

"Plasmid" => false

"Repeats" => [1, 5, 25]
"Volume" => 50

"PrimerConc" => 10

"DNAConc" => 15


passed parameters: PCRparams

6 rows × 4 columns

	Temp	Time	TimeScale	Cycles
	Int64	Any	String	Int64
1	95	1	minute	1
2	95	15	second	30
3	52	15	second	30
4	72	18	second	30
5	72	2	minute	1
6	8	00	-	1

3 rows × 6 columns

	TaqMix	Primers	Template	Water	Units	Repeats
	Int64	Int64	Int64	Int64	String	Int64
1	25	1	2	22	μ1	1
2	125	5	10	110	μ1	5
3	625	25	50	550	μ1	25

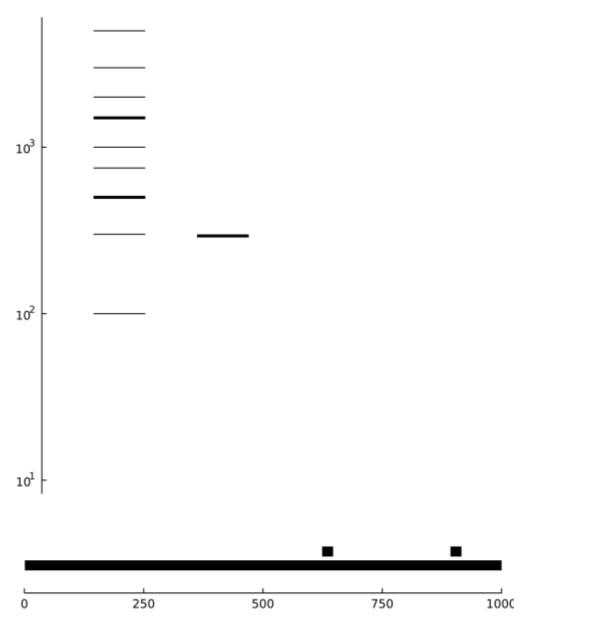

Published from PCRprep.jmd using Weave.jl v0.10.10 on 2022-04-27.